

Operator's Manual

Combination Blanket/Fluid Warmer

EC1350BL EC1850BL





MN-39940

REV.04 5/21

EN





Manufacturer's Information

Copyright

© Copyright 5/21 by Enthermics Medical Systems.

All rights reserved.

This manual or any portion thereof may not be reproduced or used in any manner whatsoever without the express written permission of Enthermics

Medical Systems.

Trademarks

All trademarks referenced in this documentation are the property of their respective owners.

Manufacturer



Enthermics Medical Systems

An ISO 13485:2016 certified company

W164 N9221 Water Street

Menomonee Falls, WI 53051, USA

Original instructions

The content in this manual is written in American English.





Authorized Representative:

MDSS GmbH Schiffgraben 41 30175 Hannover Germany



FOREWORL

Thank you for your Purchase!

This warmer has been thoroughly tested and inspected to ensure only the highest quality is provided. We supply the most durable, convenient, efficient and safe warming equipment on the market. All warmers are manufactured and fully inspected in the USA with a commitment to quality.

Register Your Warmer

Register

Register your Enthermics appliance online. Registering your appliance ensures prompt service in the event of a warranty claim. You will also receive direct notifications of software updates and additional product information.

Your personal information will not be shared with any other company.

www.enthermics.com/warranty-registration

Enthermics 24/7 Emergency Repair Service

Call

Call 800-558-8744 to reach our 24-hour emergency service call center for immediate access to local authorized service agencies outside standard business hours. The emergency service access is provided exclusively for Enthermics Medical Systems equipment and is available throughout the United States through Enthermics toll free number.

Availability

Emergency service access is available seven days a week, including holidays.



FOREWORD

This page intentionally left blank.



Manufacturer's Information	2
reword Thank you for your Purchase!	3
fety The Meaning of Signal Words	
bels Label Locations	11 11
ecifications Specification Information	
stallation How to Unpack the Warmer	
Preparing the Warmer for First Use	28 30 31 33 34
Maintenance Maintenance Schedule	41
What to do if a Power Interruption Occurs	46
ectrical Information	53
arranty Warranty	57 57



FOREWORD



SAFETY

The Meaning of Signal Words

This manual contains signal words where needed. These signal words must be obeyed to reduce the risk of death, personal injury, or equipment damage. The meaning of these signal words is explained below.



DANGER

Danger indicates a hazardous situation which, if not avoided, will result in serious injury or death.



WARNING

Warning indicates a hazardous situation which, if not avoided, could result in serious injury or death.



CAUTION

Caution indicates a hazardous situation which, if not avoided, could result in minor or moderate injury.

NOTICE

Notice indicates a situation which, if not avoided, could result in property damage.



NOTE: Note indicates additional information that is important to a concept or procedure.

Safety Precautions

Before you begin

Read and understand all instructions in this manual.

Electrical precautions

Obey these electrical precautions when using the warmer:

- Connect the warmer to a properly grounded outlet. Do not use the warmer if it is not properly grounded. Consult an electrician if there is any doubt that the outlet used is properly grounded.
- Keep the cord away from hot surfaces.
- Do not attempt to service the warmer or its cord and plug, when plugged in.
- Do not operate the warmer if it has a damaged cord or plug.
- Do not immerse the cord or plug in water.
- Do not let the cord hang over the edge of a table or counter.
- Do not use an extension cord.

Usage precautions

Obey these usage precautions when using the warmer:

- Only use this warmer for its intended use of warming medical solution bags, bottles, and/or blankets. Follow facility and solution manufacturer guidance regarding warming temperature for any item to be placed in the warmer.
- Do not use this warmer for warming blood or blood products.
- Do not cover or block any of the openings of this warmer.
- Do not use this warmer in a wet location.
- Only clean the warmer when the power cord is unplugged.
- Do not use corrosive chemicals when cleaning the warmer.
- Do not use the warmer cavity for storage.
- Do not remove exterior panels from the warmer or attempt repairs. The warmer has no user-serviceable internal components. Only perform routine cleaning and maintenance procedures specifically described in this manual. Inspection and servicing of internal components must only be performed by qualified service personnel.
- Only a qualified Enthermics service representative may make modifications to the warmer. Modifications to the warmer could be hazardous to users and patients.

Operator training

All personnel using the warmer must have proper operator training. Before using the warmer:

- Read and understand the operating instructions contained in all the documentation delivered with the warmer.
- Know the location and proper use of all controls.
- Keep this manual and all supplied instructions, diagrams, schematics, parts lists, notices, and labels with the warmer if the warmer is sold or moved to another location.



Operator qualifications

Only trained personnel with the following operator qualifications are permitted to use the warmer:

- Have received proper instruction on how to use the warmer.
- Are familiar with the purpose, limitations, and associated hazards of the warmer.

The warmer must not be used by:

People impaired by drugs or alcohol.

Condition of warmer

Only use the warmer when:

- All controls operate correctly.
- The warmer is installed correctly.
- The warmer is clean.
- The warmer labels are legible.

Servicing the warmer

- Obey precautions in the manual, on tags, and on labels attached to or shipped with the warmer.
- Only trained personnel are permitted to service or repair the warmer. Repairs that are not performed by a trained technician, or the use of non-factory parts, will void the warranty and relieve all liability.
- Any troubleshooting guides and components views included with this manual are for reference only and are intended for use by qualified and trained service technicians.
- To prevent serious injury, death or property damage, have the warmer inspected and serviced at least every twelve (12) months by a trained technician.
- Contact Enthermics for the authorized service partner in your area.

Incident notice

Any serious incident that has occurred in relation to the warmer should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.




SAFETY

This page intentionally left blank.




ABELS

Label Locations

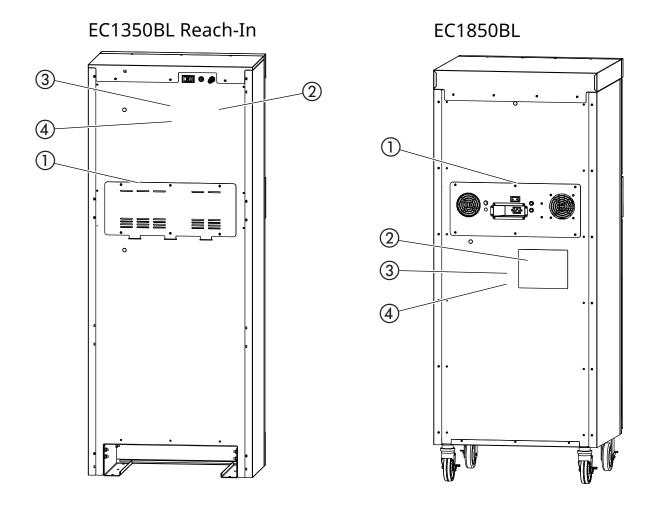

EC1350BL Pass-Through (Top View)

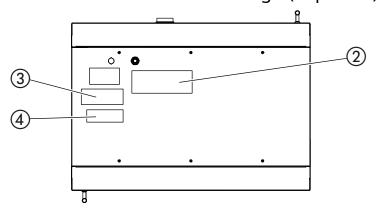

EC-PHD-004613



Security seal EC1350BL Electrical warning label LA-29320 EC1850BL 2 Electrical warning label LA-37161 **Enthermics contact information** (WarmRight. 3 ☆☆☆☆ ■ Made in ☆☆☆☆ ☆☆☆☆ ■ USA EC REP Germany Tel.: +49 511 6262 8630 Fax: +49 511 6262 8633 Quality inspected **QUALITY INSPECTED INSPECTORS INITIALS** DATE LA-36649 www.enthermics.com/quality WARNING Hot surface (label installed on the interior cavity wall)

Specification Information

Model	EC1350BL	EC1850BL	
Intended use	Bottom cavity—warming blankets TriMode cavity—warming blankets, irrigation or injection fluids		
Weight	Net: 336 lb (152 kg)	Net: 457 lb (207 kg)	
	Ship: 407 lb (185 kg)	Ship: 528 lb (239 kg)	
Storage cavity capacity	TriMode cavity: 24 1-liter bottles or 20 1-liter bags, or 9–12 blankets (3.9 ft ³)	TriMode cavity: 30 1-liter bottles or 24 1-liter bags, or 10–14 blankets (5.4 ft ³)	
сараситу	Blanket cavity: 24–30 blankets (9.4 ft ³) Blanket cavity: 30–40 blar		
Irrigation fluids (IRR): 90°F to 150°F (32°C to 66°C)			
Temperature range	Injection fluids (INJ): 90°F to 104°F (32°C to 40°C)		
	Blanket: 90°F to 160°F (32°C to 71°C)		
Clearance requirements	2" (102mm) from rear 1" (25mm) from top 1" (51mm) from sides NOTE: See topic <i>Dimension</i> <i>Drawings</i> for EC1350BL Built- In requirements.	4" (102mm) from rear 1" (25mm) from top 2" (51mm) from sides 1/2" (13mm) from bottom	

Transportation and storage conditions

- Ambient temperature range of -40°F to +159°F (-40°C to +70°C).
- Relative humidity range of 10% to 95%, non-condensing.
- Atmospheric pressure range of 7.25 psi to 15.37 psi (50kPa to 106kPa).

Operating conditions

- The warmer must acclimate to the room temperature it will be placed in—24 hours is recommended.
- The recommended environmental temperature range is 60°F to 90°F (15°C to 32°C).
- The recommended relative humidity is above 20%, non-condensing.

Standards for electrical equipment


- Medical equipment listed by Underwriters Laboratories with respect to electrical shock, fire, and mechanical hazards only, in accordance with UL 61010-1 and CAN/CSA C22.2 No. 61010-1.
- C UL US LISTED E471516
- Grounding reliability can only be achieved when the appliance is connected to an equivalent receptacle marked "Hospital Grade."
- EN61326-1; EN61000-3-2; EN61000-3-3
- IP-XO ordinary







Dimension Drawings


EC1350BL Reach-In

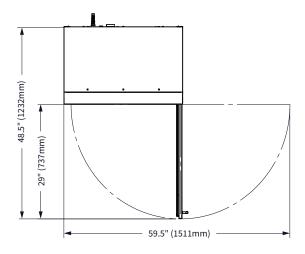


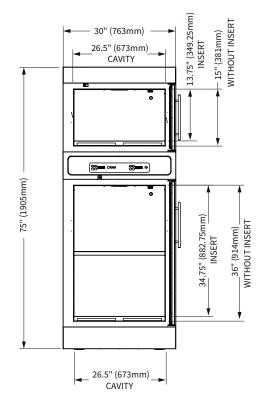


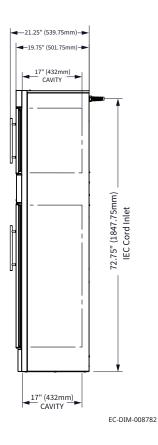

EC1350BL Pass-Through

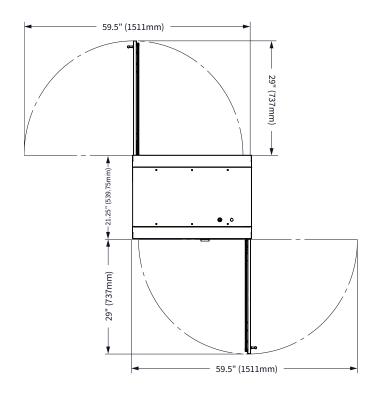


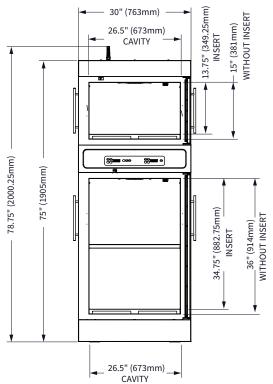


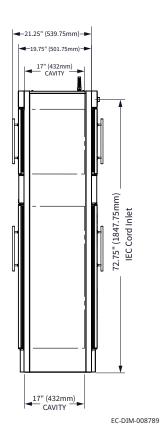




EC1350BL Built-In

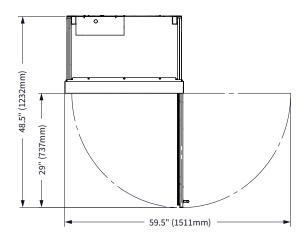

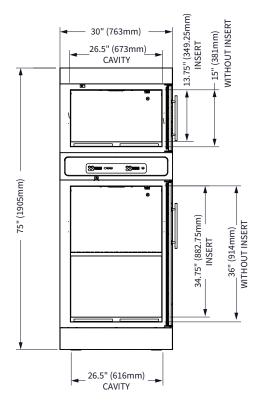

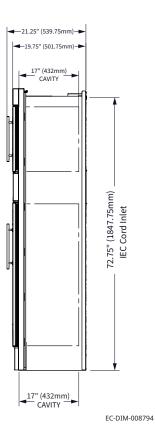

EC1350BL Built-In Cutout Specifications and Inlet Locations



CAUTION: Fire hazard.

Supply wires must be rated for a minimum wire temperature of 80°C

NOTICE

Built-in units must be hardwired to the facility mains with a 20A breaker switch.

Breaker*	Wire Size**	Wire Temperature
20A	12 AWG (1mm ²)	80°C (minimum)

^{*}Electrical connections must meet applicable codes.


Cutout Spe	cifications
Height	74.5" (1892.25mm)
Width	29.25" (743mm)
Depth	17" (432mm)
Top of unit clearance	12" (304.75mm)

	Inlet Options
1	Back of warmer
2	Top of warmer
3	Left side of warmer
_	_





^{**}IEC 60335 minimum.


Continued from previous page

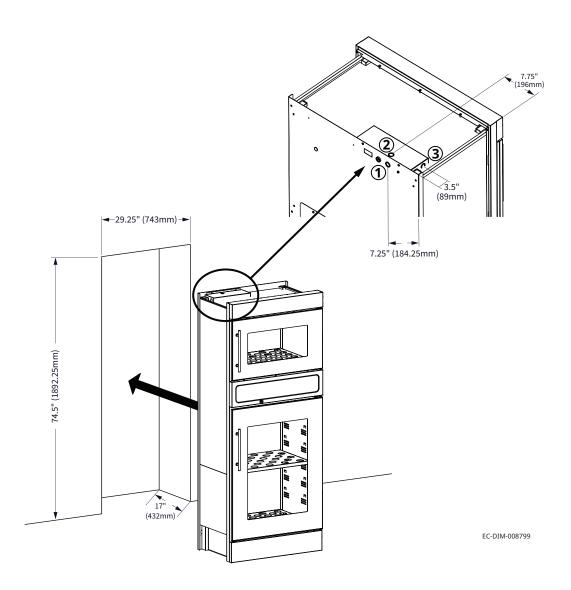


EC1850BL shown with feet

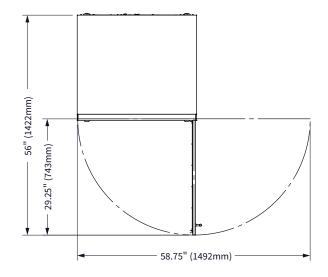


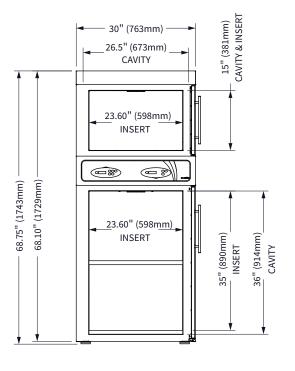


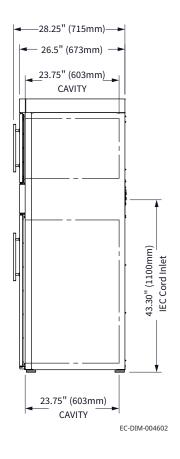

EC1850BL shown with casters

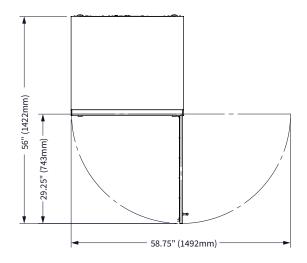

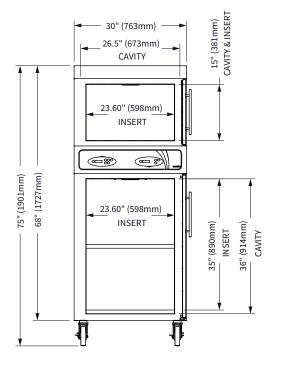

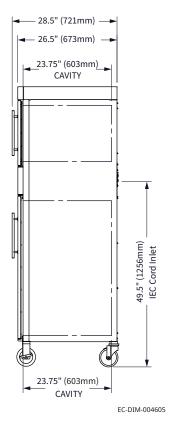

INSTALLATION

How to Unpack the Warmer



WARNING: Crushing hazard.

An unstable warmer can lead to a crushing hazard. Observe your facility's best practices for moving large equipment.

Before you begin

Make sure you have:

- Reviewed the warmer specifications.
- An appropriate lifting device and enough personnel to safely move and position the weight of the warmer.
- Cutting tools to remove the packaging.

Unpack the warmer

To unpack the warmer, do the following.

Step	Action
1.	Remove the box or crate. Save all packing materials for inspection by the carrier.
	NOTE: Examine the warmer for damage. If the warmer has been damaged, do not use the warmer until it has been inspected by an authorized service provider. Contact your carrier or customer service.
2.	Cut and remove the retaining straps and plastic wrap.
3.	Remove the warmer from the pallet using an appropriate lifting device.
4.	Remove the paperwork from the cavity.



Continued from previous page

EC1350BL only

5. Open the accessories box and install the shelf clips and shelves into the cavities.

NOTE: Set aside the toe panel(s) and screws. The toe panel(s) and screws will be used during installation after mounting the warmer.

Result

The warmer is now unpacked.

How to Install the Warmer



WARNING: Tipping hazard.

Mount the EC1350BL to the floor before use.

Before you begin

Make sure you have:

An appropriate lifting device and enough personnel to safely move and position the weight of the warmer.

Requirements

- The warmer must be installed on a level surface.
- The warmer must not be installed in any area where it may be affected by dripping water, high temperature, or any other severely adverse conditions.

Voltages

Model	v	Ph	Hz	kW
EC1350BL	120	1	50/60	1.2
	220	1	50/60	1.1
	230	1	50/60	1.1
EC1850BL	120	1	50/60	1.3
	220	1	50/60	1.2
	230	1	50/60	1.3



Continued from previous page

Position the warmer

To position the warmer, do the following.

Step Action

- 1. **Make sure** that:
 - The location where the warmer is being installed is rated to support the weight of the warmer;
 - The warmer is within five feet of the appropriate electrical outlet;
 - You follow the warmer clearance guidelines shown below. Refer to topic *Specification Information* for clearance requirements.
 - You reference the built-in specification page for cutout dimensions (EC1350BL only).
- 2. **Move** the warmer to the installation location and onto the final resting surface.
- 3. **Lock** the casters, if equipped.

The warmer is now correctly positioned.

Mount the warmer to floor (EC1350BL only)

To mount the warmer to the floor, do the following.



WARNING: Tipping hazard.

Mount the EC1350BL to the floor before use.

Non-seismic mounting hardware included in kit.




NOTE: Reference OPM-0562-19 for California seismic conformance.

4. **Bolt** the warmer to the floor.



Continued from previous page

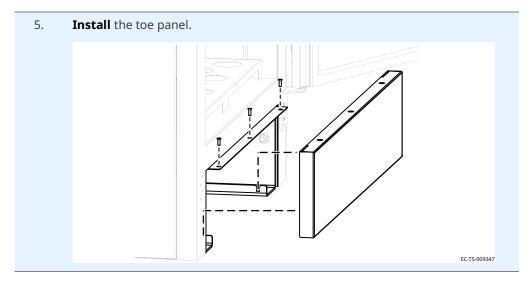

Connect power

To connect electric power to the warmer, do the following.

6. **Connect** the plug to the electrical outlet.



NOTE: For built-in units, contact a licensed contractor or electrician for connection to facility power.

Result

The warmer is now installed and ready to be used.



INSTALLATION

This page intentionally left blank.



OPERATION

Preparing the Warmer for First Use

Before you begin



CAUTION: Burn hazard.

Allow the warmer to cool before cleaning.

NOTICE

Do not use:

- abrasive cleaning compounds.
- chloride based cleaners.
- commercial or household cleaners containing ammonia.

Approved cleaning agents

- 70% or greater isopropyl alcohol
- PDI Sani-Cloth® HB

Procedure

To prepare the warmer for first use, do the following.

Step	Action
1.	Make sure that the warmer is turned off and cool.
2.	Wipe the outside of the warmer:
	with a stainless steel cleaner, if stainless steel panels.
	with an approved cleaning agent, if painted panels.
3.	Clean the interior of the warmer with a damp cloth or approved cleaning agent.
4.	Dry the interior and exterior of the warmer with a clean, lint-free cloth.
	Leave the door open until the interior of the warmer completely dries.
5.	Clean each side of the window pane with glass cleaner or distilled vinegar.

Result

The warmer is now ready for use.




How to Turn On and Turn Off the Warmer

Before you begin

The warmer must be connected to electric power.

Turning on the warmer

To turn on the warmer, do the following.

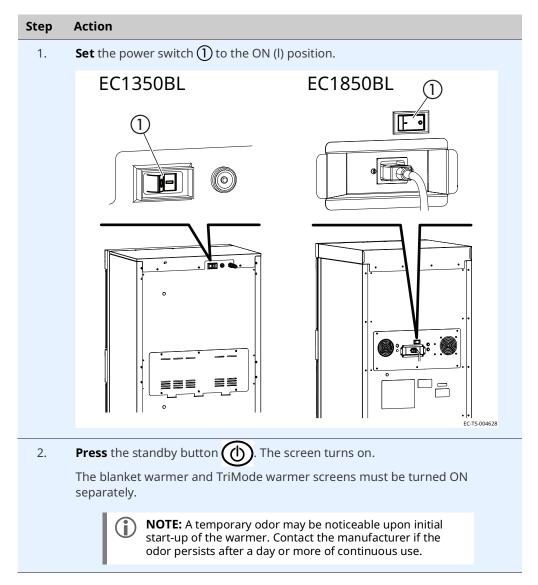

The warmer is now on.

Continued from previous page

Turning off the warmer

To turn off the warmer, do the following.

3. **Press and hold** the standby button until the screen turns off.

The blanket warmer and TriMode warmer screens must be turned OFF separately.

The warmer is now off.



How to Operate the Blanket Warmer (Bottom Cavity)

Before you begin

The warmer must be connected to electric power and turned on.



WARNING: Personal injury hazard.

Do not operate the warmer in the presence of a flammable anesthetic mixture with air, oxygen, or nitrous oxide; in oxygen-enriched environments; or in any other potentially explosive environment.



WARNING: Personal injury hazard.

Make sure the blanket support assembly and shelves are installed to prevent the blankets from being scorched or discolored.

NOTICE

Do not overload the cavity.

Allow 1" (25mm) between the blankets and the interior walls.

Procedure

To operate the warmer, do the following.

Step **Action**

1. **Set** the temperature using the arrow buttons





NOTE: The temperature set-point range is 90°F – 160°F (32°C - 71°C).

2. **Open** the door and load the dry, cotton blankets into the warmer.

Close the door.



Do not warm items containing plastic, rubber, metal snaps, studs, hooks, etc.

During the warming process

3. **Press** the temperature recall button temperature.



to view the measured cavity

The measured cavity temperature will display for five seconds. Then, the set-point temperature will display.

4. **Rotate** the blankets daily from the bottom to the top of the load to ensure equal heating and to prevent discoloration.

Result

The blankets are now warming.

How to Operate the TriMode Warmer (Top Cavity)

Before you begin

The warmer must be connected to electric power and turned on.



WARNING: Personal injury hazard.

Do not operate the warmer in the presence of a flammable anesthetic mixture with air, oxygen, or nitrous oxide; in oxygen-enriched environments; or in any other potentially explosive environment.



WARNING: Personal injury hazard.

Verify the fluid temperature prior to using the fluid. Refer to the fluid manufacturer's label for recommended warming procedures.

Do not use any fluids that are warmed above the suggested temperature.

NOTICE

Do not overload the cavity.

Refer to topic *Specification Information* for the storage cavity capacity.

Procedure

To operate the warmer, do the following.

Step Action

- 1. **Select** the operation mode.
 - Press the irrigation button for irrigation mode.
 - **Press** the injection button for injection mode.

If irrigation or injection are selected, the operation mode will illuminate



The blanket operation mode icon will not illuminate.

Press the blanket button



NOTE: Open the door to cool the cavity before changing from a higher temperature to a lower temperature.

for blanket mode.



PERATION

Continued from previous page

Set the temperature using the arrow buttons





NOTE: The irrigation temperature set-point range is 90°F – 150°F (32°C – 66°C).

The injection temperature set-point range is 90°F – 104°F (32°C - 40°C).

The blanket temperature set-point range is 90°F – 160°F

3. **Open** the door and load the fluid or blankets into the warmer.

Close the door.

During the warming process

Press the temperature recall button to view the measured cavity 4. temperature.



The measured cavity temperature will display for five seconds. Then, the set-point temperature will display.

Result

The fluids or blankets are now warming.



How to Change the Temperature Scale

Before you begin

The warmer must be connected to electric power and the screen turned off.

Procedure

To change the temperature scale from °F to °C and vice versa, do the following.

Step	Action
1.	Press and hold the standby button until the screen turns off, then release the button.
2.	Press and hold the temperature recall button until the temperature digits disappears and only the F or C display, then release the button.
3.	Press the up and down arrow buttons to toggle between the temperature scales.
4.	Press the standby button to turn the screen on.

Result

The temperature scale has now been changed.

How to Change the Sound Settings

Before you begin

The warmer must be connected to electric power.

Procedure

To change the sound settings, do the following.

Step	Action
1.	Press and while holding the temperature recall button , press the down arrow button . Once the sound volume setting displays, release both buttons.
2.	Press the up and down arrow buttons to adjust the volume. The volume range is 0 (mute) to 12 (maximum).
	The volume range is a (mate) to 12 (maximam).
3.	Press the standby button to turn the screen on.

Result

The sound setting has now been changed.

How to Lock and Unlock the Controller

Before you begin

The warmer must be turned on (screen is on).

Background

The controller can be locked to prevent changes being made to the temperature set-point.

Locking the controller

To lock the controller, do the following.

Step **Action**

Press and while holding the standby button (1), press the up arrow 1. button

The lock indicator





NOTE: Make sure the set-point temperature did not change during the locking process.

Unlocking the controller

To unlock the controller, do the following.

Press and while holding the standby button (1), press the down arrow 2. button



The lock indicator qoes off.



NOTE: Make sure the set-point temperature did not change during the unlocking process.

Result

The controller has now been locked or unlocked.



How to Download Temperature Data from the Warmer

Before you begin

- The warmer must be turned on (screen is on).
- Do not remove the USB drive during the download process.
- You will need a SanDisk Cruzer Glide USB drive.



Only use a SanDisk Cruzer Glide USB drive larger than 4 GB. Using other USB drives may cause erratic operation or loss of data.

Background

The warmer automatically downloads only the data that has not been previously downloaded. The warmer will store multiple years of data. The download time for one week of data should be under five seconds.

The warmer records various performance parameters during the following conditions:

- automatically every five minutes while the power switch is ON.
- each time the temperature set-point is changed.
- each time the TriMode cavity door is opened.
- each time the fluid operation mode is changed.

Procedure

To download temperature data from the warmer, do the following.

Step	Action
1.	Open the blanket warmer door.




Continued from previous page

2. **Plug** the USB drive into the port ①. The warmer downloads the temperature data.

NOTICE

Do not remove the USB drive during the download process.

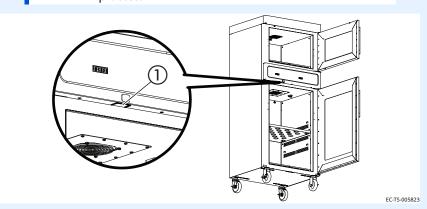

The warmer goes through the download process:

- The LED screen located above the USB port flashes for a few seconds to show a data transfer is in progress.
- The LED screen stops flashing.
- The download is complete.
- 3. **Remove** the USB drive from the port.
- 4. **Close** the blanket warmer door.

Result

The data has now been downloaded.

Viewing the temperature data

To view the temperature data, do the following.

- 5. **Plug** the USB drive into a computer.
- 6. **Double-click** on the file.



NOTE: Data is stored on the USB drive in a file named, datalog_MMDDYY_HHMMSS.csv.

7. **Sort** the data as necessary.

Temperature data will be recorded in the same temperature units that are displayed on the LED screen.



NOTE: Errors recorded by the warmer are intended for use by service technicians only.



Continued from previous page

Viewing data in Microsoft® Excel

The table below is an example of temperature data from the warmer.

Time Stamp	Chamber	Mode	Power	Door	Fan	Temperature	Set-point	Errors
12/19/2018 14:48	bottom		off	closed	off	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:48	top		off	closed	off	71.2	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:53	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:53	top	INJ	on	closed	on	70.4	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:58	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 14:58	top	INJ	on	closed	on	70.4	104	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 15:03	bottom		on	closed	on	71.7	130	Diagnostics.SensorSho rtErrors=0x0c
12/19/2018 15:03	top	INJ	on	closed	on	70.3	104	Diagnostics.SensorSho rtErrors=0x0c

Result

The temperature data has now been viewed.



OPERATION

This page intentionally left blank.



MAINTENANCE

Maintenance Schedule



NOTE: Do not remove exterior panels from the warmer or attempt repairs. The warmer has no user-serviceable internal components. Only perform routine cleaning and maintenance procedures specifically described in this manual.

Daily

For daily maintenance, do the following.

Check:

- the air vents in the airflow insert panels are not obstructed (if applicable).
- all fan guards are clear and not obstructed (if applicable).
- the number of bottles/bags, or blankets as applicable do not exceed the maximum capacity per shelf or basket.

Clean:

any spills with a clean, lint free cloth. See topic How to Clean the Warmer for the list of approved cleaners.

Monthly

For monthly maintenance, do the following.

Check:

- the door gasket for tears and holes. Make sure that it is firmly attached to the door. Check the seal when the door is closed.
- the guards around the air temperature sensor are in place and fully secure to the warmer.
- the hardware securing the warmer(s) to the wall, if applicable.

Clean:

- the outside of the warmer with stainless steel cleaner. See topic *How to Clean the* Warmer for the list of approved cleaners.
- vacuum (if applicable):
 - fan openings
 - fan sail switch
 - vent openings



MAINTENANCE

Continued from previous page

Yearly

For yearly maintenance, do the following.

Check:

- the set-point temperature compared to the actual temperature displayed.
 - Check the cavity air temperature with a thermocouple placed 1" (25mm) from the cavity sensor. Do not allow the thermocouple to touch any surface.
 Monitor the temperature for approximately one hour in an empty cavity.
- the condition of the plug and cord and replace if damaged.
- the controller screen for excessive wear. Make sure the controller screen displays and operates properly.

Clean:

• the shelves and interior of the warmer. See topic *How to Clean the Warmer* for the list of approved cleaners.



How to Clean the Warmer

Before you begin



WARNING: Electric shock hazard.

Disconnect the warmer from electric power before cleaning.



WARNING: Crushing hazard.

An unstable warmer can lead to a crushing hazard. Observe your facility's best practices for moving large equipment.



CAUTION: Burn hazard.

Allow the warmer to cool before cleaning.

NOTICE

Do not use:

- abrasive cleaning compounds.
- chloride based cleaners.
- commercial or household cleaners containing ammonia.

Approved cleaning agents

- 70% or greater isopropyl alcohol
- PDI Sani-Cloth® HB

Monthly cleaning procedure

To clean the warmer monthly, do the following.

Step	Action
1.	Make sure that the warmer is disconnected from electric power and cool.
2.	Wipe the outside of the warmer with a stainless steel cleaner.
3.	Vacuum lint from fan openings, fan sail switch, and vent openings.
4.	Clean the interior of the warmer with a damp cloth or approved cleaning agent.
5.	Dry the interior and exterior of the warmer with a clean, lint-free cloth.
	Leave the door open until the interior of the warmer completely dries.
6.	Clean each side of the window pane with glass cleaner or distilled vinegar.




Continued from previous page

Yearly cleaning procedure


To clean the warmer yearly, do the following.

Step Action

- 1. **Make sure** that the warmer is disconnected from electric power and cool.
- 2. **Remove** the shelves (1) from the interior of the warmer.

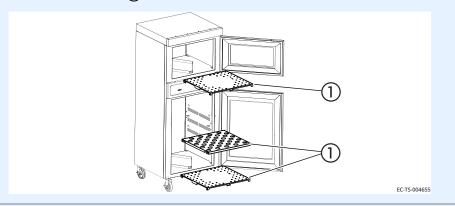

3. **Remove** the side panels ② from the interior of the blanket warmer cavity.

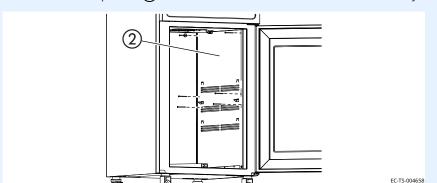

- 4. **Clean** the shelves and interior with a damp cloth or approved cleaning agent.
- Dry the interior of the warmer with a clean, lint-free cloth.
 Leave the door open until the interior of the warmer completely dries.
- 6. **Clean** each side of the window pane with glass cleaner or distilled vinegar.
- 7. **Re-install** the side panels and shelves.



WARNING: Personal injury hazard.

The blankets or fluids may overheat if the side panels and shelves are not in place. Make sure to re-install the side panels and shelves before operating the warmer.

Result

The warmer is now clean.

Warmer Disposal / Decommissioning

This product and its accessories must be disposed of according to local laws and regulations. Do not dispose of this product as unsorted municipal waste.



MAINTENANCE

This page intentionally left blank.



What to do if a Power Interruption Occurs

Background

You may need to reset the warmer in the case that a power interruption occurs. When the power is restored, the ON/OFF status indicator LED decimal flashes

The controller stores all settings and continues operation using these settings in the event of a power interruption.

Procedure

To continue operation of the warmer, do the following.

Step	Action
1.	Press the standby button . The status indicator LED decimal goes out.
	If an error code displays, refer to topic <i>Error Codes</i> .
2.	Make sure the fluids or blankets are at the correct temperature before using.
	NOTE: Inspect the fluids for their expiration date before use.

Result

Resume operation of the warmer.



What to do if the Alarm Indicator Light Flashes

Background The alarm indicator light of flashes and the controller sounds an alarm when the warmer malfunctions.

Procedure If the alarm indicator light flashes, do the following.

Step	Action
1.	Press the standby button to acknowledge and mute the alarm.
	The alarm indicator light goes off.
	Refer to the troubleshooting section for the error code.
2.	Make sure the fluids or blankets are the correct temperature before using.

Result Resume operation of the warmer.

Error Codes

Background

This section is provided for the assistance of qualified and trained service technicians only and is not intended for use by untrained or unauthorized service personnel. Failure to observe this precaution may void the warranty.



NOTE: If the warmer is not operating properly, check the following before calling an authorized service agent:

Verify that the power to the warmer is on.

If applicable, ensure the female end of plug is securely seated in the warmer and that the male end of plug is in an appropriate, functioning outlet.

If applicable, examine the fuses. Replace the fuses, refer to topic *How to Replace a Fuse*.

If applicable, examine the high limit manual reset button. If the high limit manual reset button is tripped, reset the warmer, refer to topic *How to Manually Reset the Warmer*.



NOTE: All non-critical codes can be cleared using the standby button. Critical errors (marked with a *) can only be cleared by setting the power switch at the rear of the warmer to the off (O) position and allowing the warmer to cool.

NOTICE

Do not attempt to repair or service the warmer beyond this point. Contact the manufacturer for the nearest authorized service agent. Repairs made by any other service agent without prior authorization by the manufacturer will void the warranty.

Code	Refers to	Action required
Display flashes set point	Cavity temperature higher than set point	Cavity temperature is higher than the set point temperature. Allow cavity to cool to set point temperature.
E-10 ES10 ES20	Cavity sensor Sensor 1 Sensor 2	Sensor is shorted. Software disengages heating pads. Acknowledge error by pressing the standby button. If error persists, a qualified service technician should test the sensor.
ES30 ES40 ES50 ES60 ES70	Sensor 3 Sensor 4 Sensor 5 Sensor 6	■Test the sensor. Detach the sensor from the warmer. Use an Ohm meter to measure the resistance of the sensor. Check the sensor at 25°C (77°F). If the reading is 10 KOhm ±1.5 KOhm, replace the display. If the reading is ±2 KOhm, replace the sensor.
ES70 Selisor 7		Check the wires for integrity. Inspect the connections at the control and terminal block to ensure proper and secure connections. If necessary, re-secure the faulty connections.
		■ Contact service if error persists.



TROUBLESHOOTING

Code	Refers to	Action required		
E-11 ES11 ES21	Cavity sensor Pad sensor 1 Pad sensor 2	Sensor is open. Software disengages heating pads. Acknowledge error by pressing the standby button. If error persists, a qualified service technician should test the sensor.		
ES31 ES41 ES51 ES61 ES71	Pad sensor 3 Pad sensor 4 Pad sensor 5 Pad sensor 6 Pad sensor 7	■ Test the sensor. Disconnect the sensor from the warmer control PCB and use an Ohm meter to measure the resistance of the sensor. Sensor PR-37140 should measure 10 KOhm ± 1.5 KOhm at 25°C (77°F). Sensor SN-33541 should measure 100 Ohm ± 10 Ohm at 25°C (77°F). Replace sensor if measurement is outside of tolerance.		
		■ Check the wires for integrity. Check for proper and secure connections at the control and terminal block. If necessary, re-secure the faulty connections.		
		■Contact service if error persists.		
E-30	Under temperature	The cavity temperature is lower than the set temperature for 90 minutes or longer.		
		■ Make sure the door is closed.		
		If the cavity is overloaded, redistribute the inventory. Do not exceed the height of the insert.		
		■ Test the sensor. Disconnect the sensor from the warmer control PCB and use an Ohm meter to measure the resistance of the sensor. Sensor PR-37140 should measure 10 KOhm ± 1.5 KOhm at 25°C (77°F). Sensor SN-33541 should measure 100 Ohm ± 10 Ohm at 25°C (77°F). Replace sensor if measurement is outside of tolerance.		
*E-31	Cavity sensor	■The sensor reading is above the temperature set-point. Blanket warming mode triggers at 15° over set-point. Fluid warming mode triggers at 5° over set-point.		
		■The difference between the room temperature and the fluid set-point temperature must be greater than 11°C (20°F).		
		■ Contact service if error persists.		
P131	Pad sensor 1	Heater pad over-temp error.		
P231 P331	Pad sensor 2 Pad sensor 3 Pad sensor 4 Pad sensor 5	■Software disengages heating pads.		
P431 P531		Acknowledge error by pressing the standby button.		
P631	Pad sensor 6	■ Allow the warmer to cool.		
P731	Pad sensor 7	■ Contact service if the error persists.		
*E-33	Cavity sensor	■ Sensor reading is above maximum allowable temperature set-point and over temp value. Blanket warmers trigger at 82°C (180°F). Fluid warmers trigger at 71°C (160°F).		
		■ Contact service if error persists.		
*E-50	Analog to Digital Converter Error	Remove inventory, discard if necessary, and allow the warmer to cool down.		
		■ If error persists after cool down and reset, the control assembly should be replaced by a qualified service technician. Contact service.		



Code	Refers to	Action required		
E-60	Real Time Clock Checksum Error (Blanket warmers only)	Real time clock rechargeable battery backup has discharged. Plug the warmer into the outlet for 30 minutes.		
*E-61	Real Time Clock (Blanket warmers only)	Real time clock not responding. Contact service if error persists.		
E-62	Real Time Clock (Blanket warmers only)	Timer overlay is present, but no real time clock is detected. Contact service.		
*E-70	Pad Count Error	More heater pads detected than set for. Hold the standby button for 12 seconds until display shows "PAd#" (# = number of pads selected [3-7]). Press up or down arrow to adjust to correct number of pads. [Blanket warmers: 1.5ft ³ , 2.5ft ³ = 3 pads, 3.5ft ³ , 4.0ft ³ , 7.5ft ³ = 4 pads.		
		Fluid warmers: 2.5ft ³ & 4.0ft ³ = 3 pads].		
*E-71	Personality Error	Contact service.		
E-80	EEPROM Error	EEPROM not responding. Contact service if error persists.		
*E-81	Calibration not locked	Contact service.		
*E-83	EEPROM Error	Contact service for help resetting the control.		
E-87	EEPROM Error	Stored offsets corrupted. Offsets reset to 0.		
		■Control may need to be re-calibrated.		
		■ Possible bad EEPROM.		
		■ Contact service if error persists.		
E-90	Button stuck	A button has been held down for >60 seconds. Adjust control. Error will reset when the problem has been resolved.		
E-95	Factory Test pin short detected.	Ensure that debris is not causing a short between the test pins. If the pins are clean, replace the control.		
*E-98	Temperature Delta Error	Temperature of the cavity sensors 1 and 2 differ by more than 3.3°C (6°F). Remove product and allow the warmer to cool down.		
		■ Verify that the product sensor is clean and operating correctly.		
		 Set the power switch to the OFF position to clear the error code. If error persists, the cavity sensor should be replaced by a qualified service technician. Contact service. 		
E-99	Hardware Over Temp	 Inspect the connections and condition of high limit bimetal thermostat. If applicable, make sure that the compartment fan motor is operating. Contact service if error persists. 		



TROUBLESHOOTING

Code	Refers to	Action required
*EFAn	Fan or Fan Sensor Failure	 Check to make sure the fan sensor wires did not disconnect from the fan or the control board. If the error persists after checking the wires, replace the fan.
Flashing yellow iv fluid temperature status indicator	Flashing yellow LED is illuminated within iv storage shelf	 Set the power switch to the OFF position. Remove any iv bags from inside the warmer. Set the power switch to the ON position. Visually check to see if flashing yellow LED remains illuminated. Contact service if error persists.

^{*}All non-critical codes can be cleared using the standby button. Critical errors (marked with a *) can only be cleared by setting the power switch at the rear of the warmer to the OFF position and allowing the warmer to cool.



How to Replace a Fuse (EC1850BL only)




WARNING: Fire and electrical shock hazard.

Use only UL listed 20A, 250V fast acting fuses, $6.3 \text{mm} \times 32 \text{mm}$ (F1, F2). Access should be made by qualified service technicians only.

Procedure

To replace a fuse, do the following.

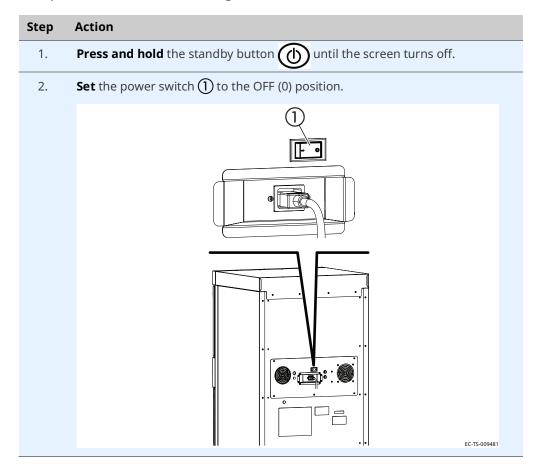




TROUBLESHOOTING

Continued from previous page

3. **Unplug** the power cord ② from the electric power source and the power inlet ③ on the warmer.

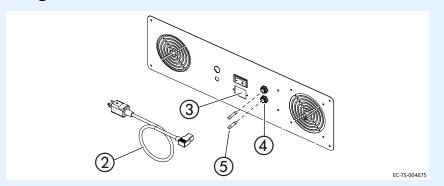

- 4. **Unscrew** the cap on the fuse holder **(4)**.
- 5. **Remove** the fuse **(5)** from the fuse holder.
- 6. **Replace** the fuse with a new fuse.
- 7. **Re-install** the cap back onto the fuse holder.

Result

The fuse has now been replaced.

Guidance and Manufacturer's Declaration

The warmer requires special precautions regarding EMC (Electromagnetic Compatibility) and needs to be installed and put into service according to the EMC information provided in the accompanying documents.

Portable and mobile RF communications equipment can affect medical electrical equipment.

A risk of increased emissions or decreased immunity may result if the power cord is altered or a manufacturer supplied power cord is not used.

The warmer should not be used adjacent to or stacked with other equipment.

The essential performance of the warmer is to not exceed an internal temperature of 160°F (71°C) for blanket warmers or 150°F (66°C) for fluid warmers.

Electromagnetic emissions

The warmer is intended for use in the electromagnetic environment specified below.

Emission test	Compliance	Electromagnetic environment - guidance			
RF emissions; CISPR 11	Group 1	The warmer uses RF energy only for internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.			
RF emissions; CISPR 11	Class B	The warmer is suitable for use in all establishments, including domestic establishments and those directly connected to			
Harmonic N/A emissions; IEC 61000-3-2		the public low-voltage power supply network that supplies buildings used for domestic purposes.			
Voltage N/A fluctuations/Flicker emissions; IEC 61000-3-3					



ELECTRICAL INFORMATION

Continued from previous page

Electromagnetic immunity

The warmer is intended for use in the electromagnetic environment specified below.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electromagnetic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines; ±1 kV for input/output lines	+2 kV for power supply lines	Main power quality should be that of a typical commercial or hospital environment. The warmer does not have any input/output lines.
Surge IEC 61000-4-5	±1 kV differential mode; ±2 kV common mode	±1 kV differential mode; ±2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5% UT (>95% dip in UT) for 0.5 cycle 40% UT (60% dip in UT) for 1 cycle 70% UT (30% dip in UT) for 25 cycles <5% UT (>95% dip in UT) for 250 sec	<5% UT (>95% dip in UT) for 0.5 cycle 40% UT (60% dip in UT) for 1 cycle 70% UT (30% dip in UT) for 25 cycles <5% UT (>95% dip in UT) for 250 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the warmer requires continued operation during power mains interruptions, it is recommended that the warmer be powered from an uninterrupted power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: UT is the a.c. mains voltage prior to application of the test level.



Continued from previous page

Electromagnetic emissions

The warmer is intended for use in the electromagnetic environment specified below.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance	
Conducted RF IEC 61000-4-6	3 V/m 150 kHz to 80 MHz	3 V/m	Portable and mobile RF communications equipment should be used no closer to any	
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2.7 GHz		part of the warmer, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.	
			Recommended separation distance	
			d = [3.5/3] √P	
			d = [3.5/3] √P 80 MHz to 800 MHz	
			d = [7/3] √P 800 MHz to 2.5 GHz	
			where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).	
			Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ¹ , should be less than the compliance level in each frequency range ² .	
			Interference may occur in the vicinity of equipment marked with the following symbol:	

NOTE: At 80 MHz and 800 MHz, the higher frequency range applies. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.



^{1.} Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the warmer is used exceeds the applicable RF compliance level above, the warmer should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the warmer.

^{2.} Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [VI] V/m.

ELECTRICAL INFORMATION

Continued from previous page

Electromagnetic immunity distance

The warmer is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The user of the warmer can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the warmer as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m			
	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2.5 GHz	
	$d = \left[\frac{3.5}{3}\right] \sqrt{P}$	$d = \left[\frac{3.5}{3}\right] \sqrt{P}$	$d = \left[\frac{7}{3}\right] \sqrt{P}$	
0.01	0.117	0.117	0.233	
0.1	0.369	0.369	0.738	
1	1.167	1.167	2.333	
10	3.689	3.389	7.379	
100	11.667	11.667	23.333	

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

WARRANTY

Warranty

Introduction

Enthermics Medical Systems warrants to the original purchaser only, that any original part found to be defective in material or workmanship will be replaced with a new or rebuilt part at Enthermics option, subject to provisions hereinafter stated.

Warranty Period

The warranty period is as follows:

- For warming cabinets shipped to the United States or Canada,
 - The labor warranty remains in effect for one (1) year from installation or fifteen (15) months from the shipping date, whichever comes first. Enthermics will bear normal labor charges performed during the standard business hours, excluding overtime, holiday rates or any additional fees.
 - The original parts warranty remains in effect five (5) years from installation of appliance or sixty-three (63) months from the shipping date, whichever comes first.
- For warming cabinets shipped outside of the United States or Canada,
 - The original parts warranty is one (1) year from the date of installation of appliance or fifteen (15) months from the shipping date, whichever comes first.
- To be valid, a warranty claim must be asserted during the applicable warranty period. This warranty is not transferable.

Exclusions

This warranty does not apply to:

- Calibration.
- Equipment damage caused by accident, shipping, improper installation or alteration.
- Equipment used under conditions of abuse, misuse, carelessness or abnormal conditions, including but not limited to, equipment subjected to harsh or inappropriate chemicals, including but not limited to, compounds containing chloride or quaternary salts, poor water quality, or equipment with missing or altered serial numbers.
- Any losses or damage resulting from malfunction, including the loss of contents or consequential or incidental damages of any kind.
- Equipment damage caused by use of any cleaning agents other than those recommended by Enthermics, including but not limited to damage due to chlorine or other harmful chemicals.
- Equipment modified in any manner from original model, substitution of parts other than factory authorized parts, unauthorized removal of any parts including legs, or unauthorized addition of any parts.



WARRANTY

Continued from previous page

■ Collateral or incidental damage as a direct result of servicing equipment built into a wall structure is not covered under warranty. It is the responsibility of the owner to bear all expense related structural repairs including, but not limited to, external electrical connections and wiring, and the removal or replacement of caulk, grout, tile, or wall covering of any kind. A service access panel for built-in equipment installations is strongly recommended.

Conclusion

This warranty is exclusive and is in lieu of all other warranties, express or implied, including the implied warranties of merchantability and fitness for a particular purpose. No person except an officer of Enthermics is authorized to modify this warranty or to incur on behalf of Enthermics any other obligation or liability in connection with Enthermics equipment.





Enthermics Medical Systems An ISO 13485:2016 certified company

An ISO 13485:2016 certified company W164 N9221 Water St Menomonee Falls, WI 53051 262-251-8356 | 800-TO-B-WARM generalinfo@enthermics.com www.enthermics.com

Made in the U.S.A