

# Comfort / Titan Series Blanket Warmers



Our blanket warmers quickly and efficiently warm blankets to precise temperatures using WarmSafe™ zone heating technology, which provides a safe even heat by adjusting temperature throughout the cabinet. Our fan-less technology is designed to provide quiet warmth for patients, allowing them to get the rest they need for recovery. And, with the absence of lint-clogged fans, these units are safe and easy to maintain.

The patient-friendly Comfort Series allows you to put aesthetically pleasing warming units right where they are needed. The lighting inside the warmers softens the patient's room, placing more emphasis on comfort and care.

The durable Titan Series, using 201&430 stainless steel, provides years of worry-free service.

**WarmSafe Zone Heating:** Heat is evenly balanced throughout the chamber to provide safe blanket warming.

**Service:** Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment:

**Glass Door:** Energy-efficient coated glass door reflects heat back into the cavity.

**Mobile:** Optional casters

**Warranty:** 3-year parts/1-year labor (US/Canada only)

**Made in the USA** with a commitment to quality

**An ISO 13485:2016 certified company**

## Standard Features

- Windowed door
- Epoxy-coated interior surfaces
- Programmable control with temperature lock feature
- Interior lighting with a comforting blue glow (Comfort Series only)
- Temperature management access point
- 350 & 750 include a heated center shelf
- Four (4) non-skid rubber feet

## Blanket Capacity



|    |                                                |
|----|------------------------------------------------|
| 6  | <b>2.5 ft³ (DC250/EC250)</b><br>4-6 blankets   |
| 12 | <b>3.5 ft³ (DC350/EC350)</b><br>10-12 blankets |
| 16 | <b>4.0 ft³ (DC400/EC400)</b><br>12-16 blankets |
| 25 | <b>7.5 ft³ (DC750/EC750)</b><br>20-25 blankets |

## Temperature



Adjustable temperature range of  
32°-71°C (90°-160°F)

## Titan Series



## Comfort Series



## Configurations (select one of each)

### Series

- ☐ Comfort (DC)    ☐ Titan (EC)

### Size

- ☐ 2.5 ft³    ☐ 3.5 ft³    ☐ 4.0 ft³  
☐ 7.5 ft³    ☐ 7.5 ft³ P/T (Titan Series only)

### Door Assembly

- ☐ Windowed door (standard)    ☐ Solid door (Titan Series only)

### Door Swing

- ☐ Right hinged (standard)    ☐ Left hinged

### Door Lock

- ☐ No lock (standard)    ☐ Manual lock

### Electrical

- ☐ 120V 1PH    ☐ 220V 1PH    ☐ 230V 1PH

## Accessories (select all that apply)

**Bumper and Caster Assembly** 3" (76mm) casters, 2 rigid, 2 swivel with Brake; overall height increases by 5" (127mm)

- ☐ DC250 or DC350 (5017546)    ☐ EC400 or EC750 (5017549)  
☐ DC400 or DC750 (5017547)    ☐ EC750/PT (5018180)  
☐ EC250 or EC350 (5017548)

**Caster Package with Plate** 3" (76mm) casters, 2 rigid, 2 swivel with Brake; overall height increases by 5" (127mm)

- ☐ DC250 or DC350 (5017542)    ☐ EC400 or EC750 (5017545)  
☐ DC400 or DC750 (5017543)    ☐ EC750/PT (5018124)  
☐ EC250 or EC350 (5017544)

### Anchorage (4 required)

- ☐ Anti-tip bracket (1014403)




### Stacking Kits

It is recommended that fluid warming units be placed at or near the bottom of the stacked configuration whenever possible for maximum stability - not available with casters




- ☐ DC250 or DC350 over DC250 or DC350 (5013021)  
☐ DC400 or DC750 over DC400 or DC750 (5013024)  
☐ EC250 or EC350 over EC250 or EC350 (5017272)  
☐ EC400 or EC750 over EC400 or EC750 (5017273)

# Comfort/Titan Series Blanket Warmers




EC/DC250

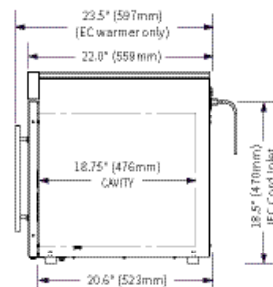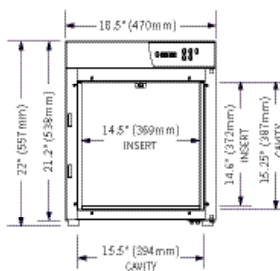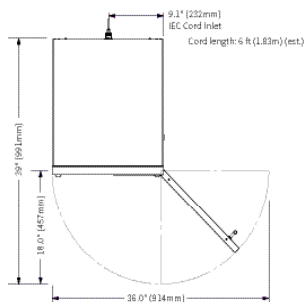

EC/DC350

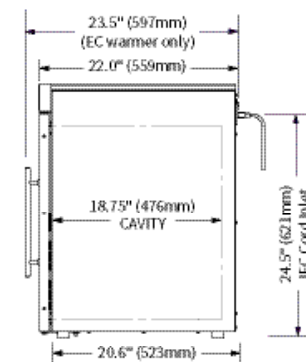
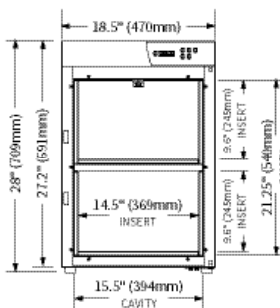
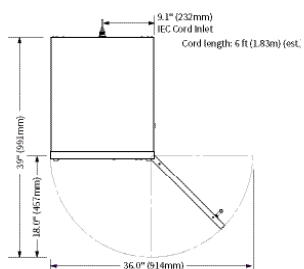


EC/DC400

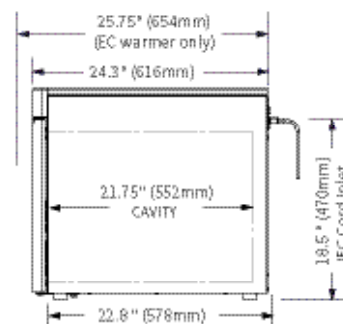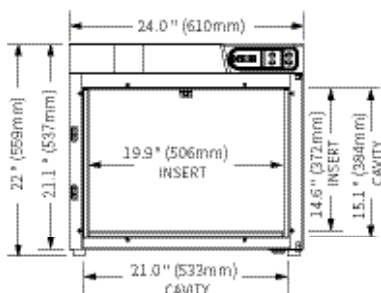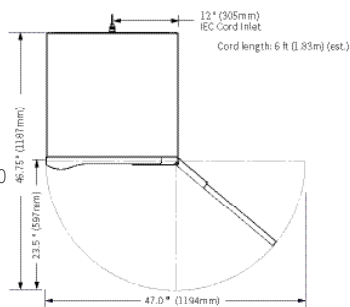

DC-004-00932

EC/DC750

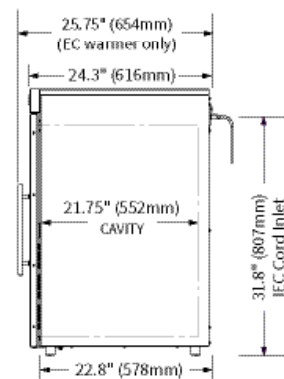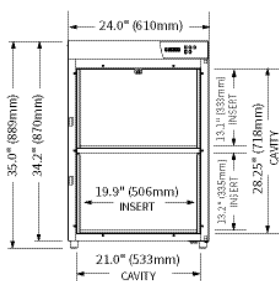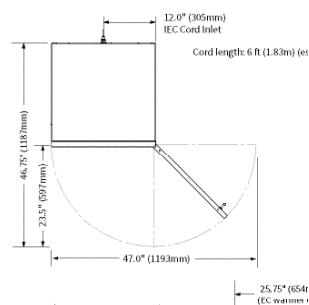

# Comfort/Titan Series Blanket Warmers

## Specification



CLEARANCE

Top: 1"  
Left: 1"  
Right: 1"  
Back: 3"



DIMENSIONS

| Model                               | Exterior (H x W x D)                       | Interior (H x W x D)                         | Net Weight      |
|-------------------------------------|--------------------------------------------|----------------------------------------------|-----------------|
| DC250                               | 22"(559mm) x 18.5"(470mm) x 22.0"(559mm)   | 15.25"(387mm) x 15.5"(394mm) x 18.75"(476mm) | 67 lbs (30 kg)  |
| EC250                               | 22"(559mm) x 18.5"(470mm) x 23.5"(597mm)   | 15.25"(387mm) x 15.5"(394mm) x 18.75"(476mm) | 67 lbs (30 kg)  |
| DC350                               | 28"(711mm) x 18.5"(470mm) x 22.0" (559mm)  | 21.25"(540mm) x 15.5"(394mm) x 18.75"(476mm) | 91 lbs (41 kg)  |
| EC350                               | 28"(711mm) x 18.5"(470mm) x 23.5" (597mm)  | 21.25"(540mm) x 15.5"(394mm) x 18.75"(476mm) | 91 lbs (41 kg)  |
| DC400                               | 22"(559mm) x 24.0"(610mm) x 24.3" (617mm)  | 15.25"(387mm) x 21.0"(533mm) x 21.75"(552mm) | 86 lbs (39 kg)  |
| EC400                               | 22"(559mm) x 24.0"(610mm) x 25.75" (654mm) | 15.25"(387mm) x 21.0"(533mm) x 21.75"(552mm) | 86 lbs (39 kg)  |
| DC750                               | 35"(889mm) x 24.0"(610mm) x 24.3" (617mm)  | 28.25"(718mm) x 21.0"(533mm) x 21.75"(552mm) | 134 lbs (61 kg) |
| EC750                               | 35"(889mm) x 24.0"(610mm) x 25.75" (654mm) | 28.25"(718mm) x 21.0"(533mm) x 21.75"(552mm) | 134 lbs (61 kg) |
| <b>Ship Dimensions (L x W x H)*</b> |                                            | <b>Ship Weight*</b>                          |                 |
| EC/DC250                            | 23" x 35" x 30"                            | 110 lbs                                      |                 |
| EC/DC350                            | 23" x 35" x 42"                            | 140 lbs                                      |                 |
| EC/DC400                            | 31" x 31" x 32"                            | 144 lbs                                      |                 |
| EC/DC750                            | 35" x 38" x 45"                            | 184 lbs                                      |                 |

\*Domestic ground shipping information. Contact factory for export weight and dimensions.



ELECTRIC

|            | V   | Ph | Hz    | A   | kW  | Plug Configuration                        |
|------------|-----|----|-------|-----|-----|-------------------------------------------|
| <b>250</b> | 120 | 1  | 60    | 5   | 0.6 | NEMA 5-15P; 15A 125V Plug; Hospital Grade |
|            | 220 | 1  | 50/60 | 2.7 | 0.6 | CEE 7/7 220-230V Plug                     |
|            | 230 | 1  | 50/60 | 2.6 | 0.6 | BS 1363 Plug; CEE 7/7 220-230V Plug       |
| <b>350</b> | 120 | 1  | 60    | 6.7 | 0.8 | NEMA 5-15P; 15A 125V Plug; Hospital Grade |
|            | 220 | 1  | 50/60 | 3.6 | 0.8 | CEE 7/7 220-230V Plug                     |
|            | 230 | 1  | 50/60 | 3.5 | 0.8 | BS 1363 Plug; CEE 7/7 220-230V Plug       |
| <b>400</b> | 120 | 1  | 60    | 5.0 | 0.6 | NEMA 5-15P; 15A 125V Plug; Hospital Grade |
|            | 220 | 1  | 50    | 2.7 | 0.6 | CEE 7/7 220-230V Plug                     |
|            | 230 | 1  | 50    | 2.6 | 0.6 | BS 1363 Plug; CEE 7/7 220-230V Plug       |
| <b>750</b> | 120 | 1  | 60    | 6.7 | 0.8 | NEMA 5-15P; 15A 125V Plug; Hospital Grade |
|            | 220 | 1  | 50/60 | 3.6 | 0.8 | CEE 7/7 220-230V Plug                     |
|            | 230 | 1  | 50/60 | 3.5 | 0.8 | BS 1363 Plug; CEE 7/7 220-230V Plug       |

\*Other international plugs are available, contact the manufacturer for more information.

\*\*Electrical connections must meet all applicable federal, state and local codes.



IPX-0

## CONTACT US

4000 S.E. Columbia Way | Vancouver, Washington 98661 | U.S.A.  
Phone: 800-822-3501 | [enthermics.com](http://enthermics.com)