

Fluid & blanket warmers

Durable, reliable temperature control.

tri mode



EC550BL-18

OVERALL SIZE
30"W x 21.12"D x 25.12"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 4.75 cubic ft. 26.5"W x 17"D x 18.25"H

- **tri mode Range:**
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 4.75 cubic ft.
- Caster (option)

tri mode



EC550BL-24

OVERALL SIZE
30"W x 28.5"D x 25.12"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 5.5 cubic ft. 26.5"W x 19.75"D x 18.25"H

- **tri mode Range:**
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 5.5 cubic ft.
- Caster (option)

tri mode



EC250L

OVERALL SIZE
18.5"W x 23.5"D x 22"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 2.5 cubic ft. 15.5"W x 18.75"D x 15.25"H

- **tri mode Range:**
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 2.5 cubic ft.
- Glide out basket
- Caster (option)

tri mode



EC350L

OVERALL SIZE
18.5"W x 23.5"D x 28"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 3.5 cubic ft. 15.5"W x 18.75"D x 21.25"H

- **tri mode Range:**
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 3.5 cubic ft.
- Glide out basket
- Caster (option)



enthermics®

by

PEDIGO®

warming cabinets

blanket | fluid | combination



EC250

OVERALL SIZE
18.5"W x 23.5"D x 22"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 2.5 cubic ft. 15.5"W x 18.75"D x 15.25"H

- **Temp Range:**
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 2.5 cubic ft.
- No fans or lint traps
- Caster (option)



EC350

OVERALL SIZE
18.5"W x 23.5"D x 28"H
COMPARTMENT SIZE
<ul style="list-style-type: none"> 3.5 cubic ft. 15.5"W x 18.75"D x 21.25"H

- **Temp Range:**
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 3.5 cubic ft.
- No fans or lint traps
- Caster (option)



EC750

OVERALL SIZE
24"W x 25.75"D x 35"H 24"W x 25.75"D x 42"H <small>(with casters)</small>
COMPARTMENT SIZE
<ul style="list-style-type: none"> 7.5 cubic ft. 21"W x 21.75"D x 28.25"H <small>* (13.18"H above and below the center shelf)</small>

- **Temp Range:**
 - Blanket: 90°F - 160°F
- **Total Capacity:**
 - 7.5 cubic ft.
- No fans or lint traps
- Caster (option)



WarmSafe™ Technology

Enthermics by Pedigo warmers feature **WarmSafe™ Zone Heat**, a method for distributing a balance of heat throughout the cabinet, to assure accurate set-point temperatures.

Fluid warmers incorporate the use of a circulation fan to evenly distribute the warmed air to achieve and maintain accurate set-point temperatures. Pedigo fluid warmers have an accuracy of +0/-2°F for injection fluids and an accuracy of +0/-3°F for irrigation fluids.



enthermics®

by

PEDIGO®

TEL: 800.822.3501 • www.enthermics.com



4000 SE Columbia Way • Vancouver, WA 98661

04/26 EBP-Trifold



TEL: 800.822.3501 • www.enthermics.com • 4000 SE Columbia Way • Vancouver, WA 98661





EC1350BL

OVERALL SIZE

- 30"W x 21.25"D x 75"H

COMPARTMENT SIZE

- 3.9 cubic ft. (top)
- 26.5"W x 17"D x 15"H
- 9.4 cubic ft. (bottom)
- 26.5"W x 17"D x 36"H

- **Temp Range:**
 - Blankets: 90°F - 160°F
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F

- **WarmSafe™ Technology**
- **Total Capacity:**
 - 13.3 cubic ft.
- Two compartments, each with individual controls
- **tri▼mode** temperature control for top compartment only
- Programmable control with auto-lock feature

EC1540

OVERALL SIZE

- 22"W x 34.5"D x 72.5"H

COMPARTMENT SIZE

- 7.7 cubic ft. (top)
- 19"W x 26.5"D x 26.75"H
- 7.7 cubic ft. (bottom)
- 19"W x 26.25"D x 26.75"H

- **Temp Range:**
 - Blankets: 90°F - 160°F
 - Irrigation: 90°F - 150°F

- **WarmSafe™ Technology**
- **Total Capacity:**
 - 15.4 cubic ft.
- Two compartments, each with individual controls
- No fans or lint traps
- Programmable control with auto-lock feature
- 5" casters

EC1850BL

OVERALL SIZE

- 30"W x 28.5"D x 75"H

COMPARTMENT SIZE

- 5.4 cu. ft. fluid
- 26.5"W x 23.75" D x 15"H
- 13.1 cu. ft. blanket
- 26.5"W x 23.75"D x 36"H

- **tri▼mode Range:**
 - Blankets: 90°F - 160°F
 - Irrigation: 90°F - 150°F
 - Injection: 90°F - 104°F

- **WarmSafe™ Technology**
- **Total Capacity:**
 - 18.5 cubic ft.
- Two compartments, each with individual controls
- **tri▼mode** temperature control for top compartment only
- 5" casters

EC2060

OVERALL SIZE

- 32"W x 29.25"D x 71"H

COMPARTMENT SIZE

- 10.3 cubic ft. (top)
- 26.5"W x 23.75"D x 28.37"H
- 10.3 cubic ft. (bottom)
- 26.5"W x 23.75"D x 28.37"H

- **Temp Range:**
 - Blankets: 90°F - 160°F

- **WarmSafe™ Technology**
- **Total Capacity:**
 - 20.6 cubic ft.
- One compartment, two separate doors
- Programmable control with auto-lock feature
- 5" casters

WarmSafe™ Technology

Zone heat warms safely, quietly and evenly.



DC250

OVERALL SIZE

- 18.5"W x 22"D x 22"H

COMPARTMENT SIZE

- 2.5 cubic ft.
- 15.5"W x 18.75"D x 15.25"H

- **Temp Range:**
 - Blankets: 90°F - 160°F
- **Total Capacity:**
 - 2.5 cubic ft.
- A warming shut-off system, separate from the electronic control, prevents overheating
- Interior blue lighting
- No fans or lint traps

DC350

OVERALL SIZE

- 18.5"W x 22"D x 28"H

COMPARTMENT SIZE

- 3.5 cubic ft.
- 15.5"W x 18.75"D x 21.25"H*
- * (9 .56"H above & below the center shelf)

- **Temp Range:**
 - Blankets: 90°F - 160°F
- **Total Capacity:**
 - 3.5 cubic ft.
- A warming shut-off system, separate from the electronic control, prevents overheating
- Interior blue lighting
- No fans or lint traps

DC750

OVERALL SIZE

- 24"W x 24.25"D x 35"H

COMPARTMENT SIZE

- 7.5 cubic ft.
- 21"W x 21.75"D x 28.25"H*
- * (13.12"H above & below the center shelf)

- **Temp Range:**
 - Blankets: 90°F - 160°F
- **Total Capacity:**
 - 7.5 cubic ft.
- A warming shut-off system, separate from the electronic control, prevents overheating
- Interior blue lighting
- No fans or lint traps

ivNow®



ivNow-1 Capacity One (1) bag, up to 3l *


ivNow-2 Capacity Two (2) bags, up to 3l *

ivNow-3 Capacity Three (3) bags, up to 3l *

The ivNow® may be placed directly on the countertop or fitted with an optional bracket for mounting on a stand or wall.

- Can provide from 25 to 78 liters of warmed fluid during a 12-hour period.
- **Easy-to-Use:** Place the IV bag in the pod to activate the unit. Fluid is automatically warmed up to 104°F (40°C) in 30 minutes or less
- **Capacity:** 1-3 pods (holds up to 3-liter bag size)
- **Configurations:** Countertop, mobile equipment pole stand, or wall.

ivNext®

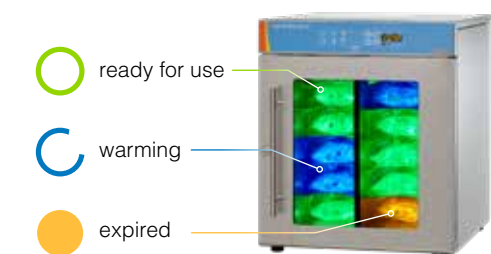


Intelligent Fluid Warming

- Know when an IV bag has reached the recommended set point temperature (adjustable injection mode)
- Color-coded lighting indicates bag readiness or expiration.
- Organized inventory allows for faster retrieval of fluids.
- **Overall Size:** 18.5"W x 23.5"D x 22"H
- **Capacity:** 2.5 cubic ft. | 10 1-liter bags

* Check with your distributor for availability if purchasing outside of the USA.

WarmSafe™ Technology

You do not need statistics or charts to tell you that patients respond well to warmed blankets and fluids. You have seen it first hand time and time again. With our patented WarmSafe™ warming technologies, our uncompromising quality construction and our variety of sizes and features, we help you maintain that kind of positive patient response – the kind that helps reduce recovery time and ensures a comfortable stay.

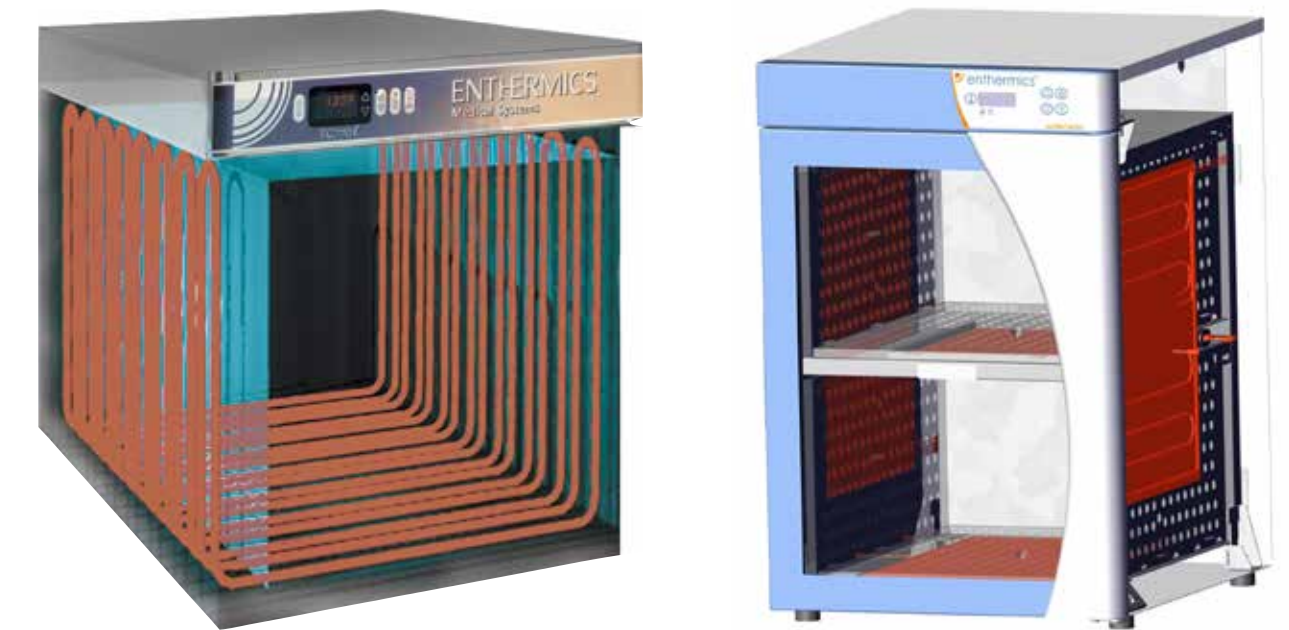

Patented warming technology distributes warmth evenly to fluids, through an electrothermal cable array in the two sides and bottom of the cabinet.

System monitors temperature in each cabinet zone and provides heat when it is needed, where it is needed. Zone heating prevents overheating and discoloration of blankets.

* available as a built-in or pass-through cabinet