

# EC1540



## Warming Cabinet

The programmable, quiet and affordable EC1540 Warming Cabinets are built for large capacity patient warming needs. The warming units come with a timer designed for energy efficiency while maintaining necessary temperature levels.

The durable, reliable stainless steel Titan Series Warming Cabinets quickly and efficiently heat blankets and irrigation fluids (temperature to be locked via digital controls to fluid manufacturer recommendation) meeting the Association of periOperative Registered Nurses (AORN) recommended temperatures and safety practices.

Our fan-less technology is designed to provide quiet warmth for patients, allowing them to get the rest they need for recovery. And, with the absence of lint-clogged fans, these units are safe and easy to maintain.

**WarmRight:** Provides the right temperature for the right product.

**Durable:** The affordable, rugged stainless steel cabinets are made to last and work just as hard as you do.

**Glass Door:** Energy-efficient coated glass door reflects heat back into the cavity.

**Mobile:** Rugged casters are designed to easily transport warmed fluids and blankets throughout the facility.

**Service:** Our worldwide support team is dedicated to providing quality equipment. We walk with you every step of the way to maximize your return on investment.

**Made in the USA** with a commitment to quality

**Warranty:** 5-year parts/1-year labor (US/Canada only)

**An ISO 13485:2016 certified company**

### Standard Features

- Windowed door
- Stainless steel exterior & interior surfaces
- Fixed shelf
- Programmable control with temperature lock feature
- Timer feature
- Temperature management access point
- Audible over-temp alarm
- Four casters - two (2) rigid and two (2) swivel with locking brake



**50** 40-50 Blankets  
**60** 1-ltr bottles

CAPACITY



TEMPERATURE

Adjustable temperature range of  
37°-82°C (98°-180°F)

**When warming irrigation fluids, lock the max temperature per fluid manufacturer recommendation**



### Configurations (select one)

#### Door Assembly

- Window door (standard)  Solid Door

#### Door Swing

- Right hinged

#### Door Lock

- No lock (standard)  Manual lock

#### Electrical

- 120V 1PH

### Accessories (select all that apply)

#### Legs

- Leg kit, 6" (152mm) (5011149)

#### Bumper

- Full-perimeter rubber bumper assembly (5015440)

#### Seismic Mounting Brackets

Special installation requirements apply

- Contact factory


#### Built in Trim Kit

- Trim kit

# EC1540 Blanket Warming Cabinet



## Specification

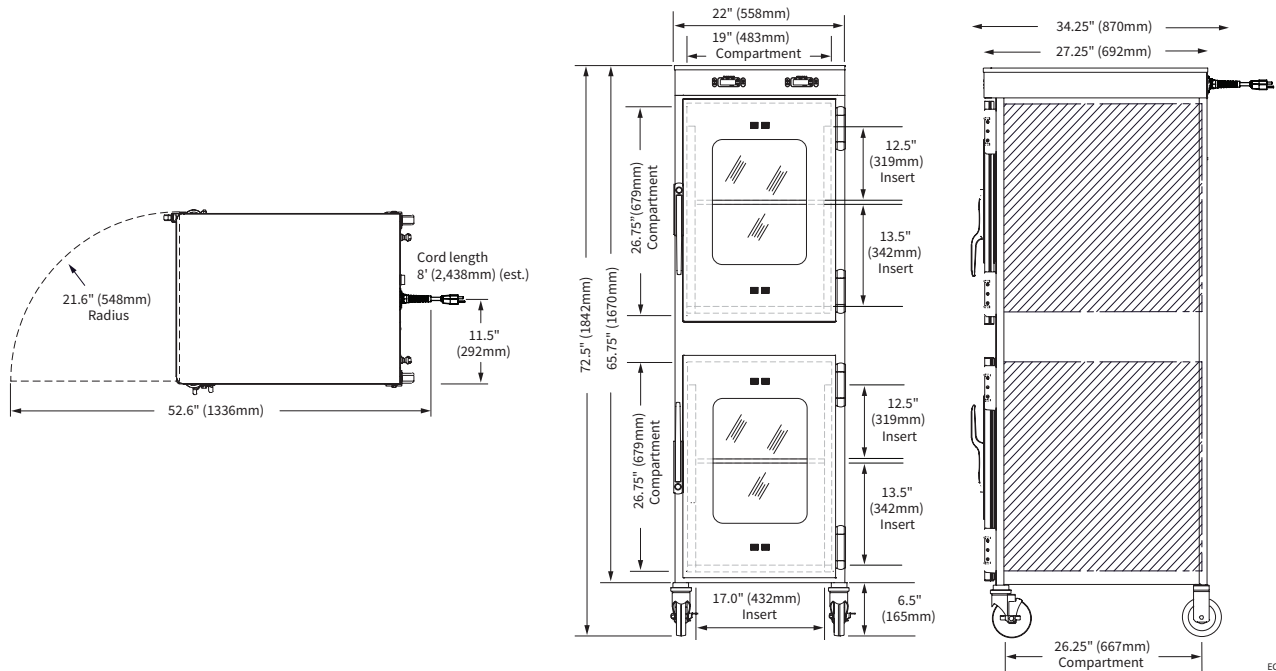

EC-DIM-005116



**Top: 1"**  
**Left: 1"**  
**Right: 1"**  
**Back: 3"**

**CLEARANCE**



Model	Exterior (H x W x D)	Interior (H x W x D)	Net
EC1540	72.5"(1842mm) x 22"(558mm) x 34.25"(870mm)	53.5"(1360mm) x 19"(483mm) x 26.25"(667mm)	374 lbs (170 kg)

**DIMENSIONS**

	Ship Dimensions (L x W x H)*	Ship Weight*
EC1540	38" x 35" x 80"	438 lbs (216 kg)

\*Domestic ground shipping information. Contact factory for export weight and dimensions.



	V	Ph	Hz	A	kW	Plug Configuration
EC1540	120	1	50/60	16	1.92	NEMA 5-20P; 20A 125V Plug; Hospital Grade

\*Electrical connections must meet all applicable federal, state and local codes.



## CONTACT US

W164 N9221 Water Street | Menomonee Falls, Wisconsin 53051 | U.S.A.  
 Phone: 262.251.8356 | 800.862.9276 | Fax: 262.251.7067 | [enthermics.com](http://enthermics.com)